Instructions For Use
Please Read Carefully

Contents

01.	Introduction	2
02.	Intended Use	2
03.	Microcurrent	3
04.	Contents	4
05.	Safety Guidelines	5
06.	Before Use	8
07.	Facial Strokes	9
08.	Using Uniwand	15
09.	Disposal	17
10.	Troubleshooting	18
11.	Warranty	19
12.	Technical Specification	20
13.	Symbols Table	21

01. Introduction

Uniwand brings the professional beauty therapy treatment to your home. Uniwand is a non-invasive facial toning tool to be used alongside microcurrent facial toning devices, such as the Uniglo, to reduce visible signs of ageing. Uniwand works by delivering low-level electrical impulses via the dual sphere treatment head to the facial muscles. The wand offers precise and flexible control, allowing users to target the specific areas of need. Microcurrent therapy should be completely pain-free. The current delivered during the stimulation is so small that in most cases, no sensation will be felt.

02. Intended Use

Uniwand is an accessory designed to be used with TensCare beauty products in the home environment to reduce signs of ageing. It is suitable for use by all who can control the device and understand the instructions.

Do not use the device for any purpose other than this intended use. This device is intended for cosmetic use on the face only.

Warning: Not suitable for use by children without adult supervision.

Important: The device should only be used with microcurrent programmes such as Uniglo FACE 2. Only use Uniwand with microcurrent parameters recommended in this manual.

Warning: If using Uniwand with the Uniglo Face and

Body Toning device, do not use FACE 1 programme or any other programme than FACE 2.

03. Microcurrent

Microcurrent stimulation delivers low current into the cells of the body. Microcurrent has been regularly used by beauty professionals to offer a non-invasive and safe face lift alternative to their clients. By stimulating the facial muscles, regular microcurrent treatments can positively impact the appearance of your skin.

Microcurrent can have preventative and corrective effects. The stimulation can help tighten skin and soften fine lines and wrinkles. The stimulation can also help muscle and tissue to be in a better state of repair to receive the body's own natural bio-electrical impulses; hence a preventative result can also be achieved.

Benefits of Microcurrent

- Encourages collagen production to plump out the tissue
- Helps to tighten and soften lines and wrinkles
- Boosts sagging muscles
- Reduces oedema: puffiness under the eyes is reduced as is dark circles under the eyes.
- Gives a youthful appearance to the skin
- Tightens and reduces open pores
- Enhances colours and texture of the skin


04. Contents

1 x Uniwand

1 x Unigel Hyaluronic Conductive Gel

1 x Instructions for Use

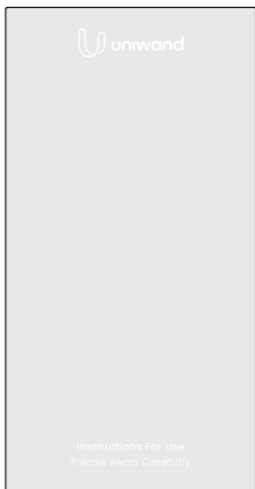05. Safety Guidelines

Contraindications

- **Do not use** if you have a pacemaker (or if you have a heart rhythm problem) or with any electronic devices. Using this unit with electronic devices may cause erroneous operation of the device. Stimulation in the direct vicinity of an implanted device may affect some models.
- **Do not use** during the first three months of pregnancy.

Warnings

- Uniwand is designed to be used with Microcurrent programmes only. Do not use with programmes that exceed the parameters detailed in the technical specification.
- Do not use when driving, operating machinery, or similar actions needing muscular control. Damaged leads, or sudden changes in contact may cause brief involuntary muscle movements.
- Do not use if you have, in the area being treated: active or suspected cancer or undiagnosed pain with a history of cancer. Stimulation directly through a confirmed or suspected malignancy should be avoided as it may stimulate growth and promote spread of cancer cells.
- Do not use while simultaneously connected to high frequency surgical equipment as it may result in burns at the site of stimulator electrodes

and possible damage to the stimulator.

- Do not use if you have metal implants in the head or neck area.
- Do not use if you have Botox injections in the application area.
- Do not use for 6-8 weeks following collagen injections.
- Do not use for 1 year following a Retin A or Ro-Accountaine treatment. This medication used in the treatment of acne can cause thinning of the skin.
- Ensure the device is switched off before inserting the Uniwand lead.

Cautions

- Uniwand is intended for single patient use only. Do not share your Uniwand. Improper treatment or cross-infection may occur.
- Caution should be used if you have a bleeding disorder as stimulation may increase blood flow to the stimulated region.
- Caution should be used if you have suspected or diagnosed epilepsy as electrical stimulation may affect seizure threshold.
- Do not permit use by children unable to understand the instructions or persons with cognitive disabilities, i.e., Alzheimer's disease or dementia.
- Keep out of reach from children. Lead wire is a long cord - risk of strangulation in infants.

- Do not use on broken skin.
- Some patients may experience skin irritation or hypersensitivity due to the electrical stimulation or electrically conductive medium (e.g., the Unigel). The irritation can usually be reduced by using an alternate conductive medium. If a skin irritation continues, discontinue treatment immediately.
- Do not immerse your device in water or place it close to excessive heat.
- Keep the device away from sunlight, as long-term exposure may affect the rubber causing it to become less elastic and crack.
- Keep the device away from lint and dust, as long-term exposure may affect the sockets or cause the battery connector to develop a bad contact.
- There are no user serviceable parts. Do not attempt to open or modify the unit. This may affect the safe operation of the unit and will invalidate the warranty.

Adverse Reactions

If you experience skin irritation or burning, discontinue use. If you experience any persistent painful sensations during or following the use of the Uniwand device, you should stop using the device and consult your healthcare professional.

06. Before Use

Uniwand has designed to be used with Uniglo's FACE 2 programme.

Important: Observe all the Contraindications, Cautions and Warnings included with the device that you are using the Uniwand with.

Programme Information

Uniglo's **FACE 2** programme has been designed to provide three phases of treatment.

Phase 1 – Circulation: This part of the treatment stimulates the epidermal and dermal blood circulation. It also has an activating effect on the lymph circulation.

Phase 2 – Superficial lift: This part of the programme lifts the superficial facial muscles.

Phase 3 – Deep facial toning: This phase gives a lift to deep facial muscles.

The programme phases are indicated on the Uniglo display by the central bar (illustration below).

The intensity levels are preset for each phase. These can be reduced by pressing the ▼ button.

First Use

Before using for the first time, test the Unigel with your skin. Apply a generous amount of gel onto any part of your body and leave for 30 minutes. After this time, clean off the gel and check the area of skin. If there is any irritation or if an allergic reaction develops within 24 hours, please discontinue use and seek advice from your healthcare professional.

Preparation

Before use ensure your skin is clean and dry. Remove facial hair as hair can reduce the conductivity between the spheres and the skin and have a negative impact on performance. Use Uniwand in front of a mirror to allow for accurate completion of the facial strokes.

07. Facial Strokes


The following pages detail microcurrent facial techniques. Each area of the face has several different strokes and starting positions.

Each facial stroke should take 5 seconds to perform.

The user should focus on one area of the face at a time.


Jaw and Cheeks

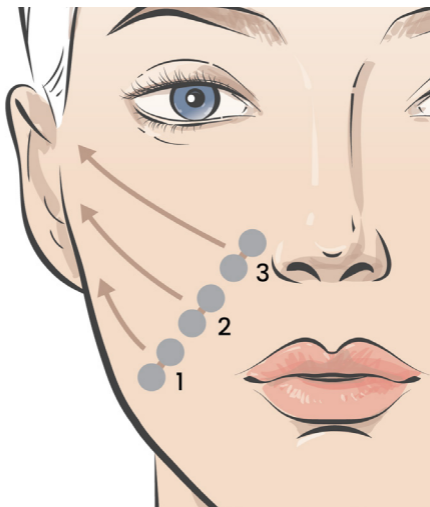- **Position 1 (Jaw)** – Position the bottom sphere on your jawline. Slowly glide the wand towards your lower ear.
- **Position 2 (Lower Cheeks)** – Position the bottom sphere under your cheekbone. Slowly glide the wand to your mid-ear.
- **Position 3 (Upper Cheeks)** – Position the bottom sphere on your cheekbone. Slowly glide the wand to your upper ear.
- Repeat for the opposite side of the face.


Forehead

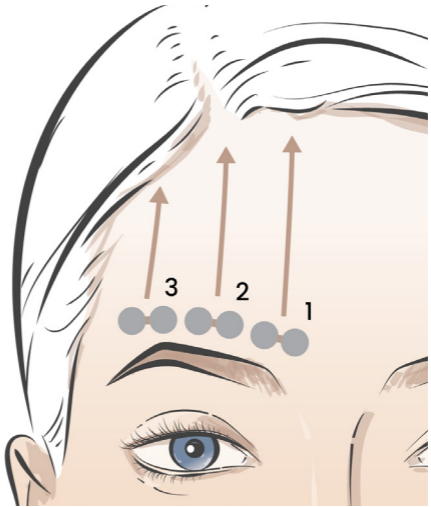- Position the wand at position 1, above the thick of your eyebrow. Slowly glide the wand toward your hairline.
- Repeat for position 2 and 3.
- Repeat for the opposite side of the face.

Eye

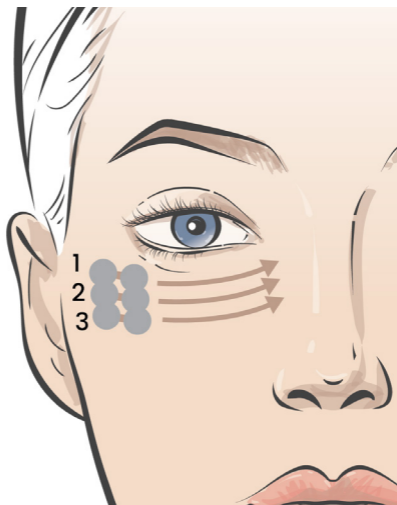
- Position the Uniwand at position 1, with the spheres positioned horizontally. Slowly glide the wand towards the inner corner of your eye, maintaining the same orientation.
- Repeat for position 2 and 3.
- Repeat for the opposite side of the face.

Lip

- Position the Uniwand above the upper lip. Push the Uniwand up towards your ear and hold the Uniwand in this position for 5 seconds.
- Repeat at position 2.
- Repeat for the opposite side of the face.

Cheek Contour

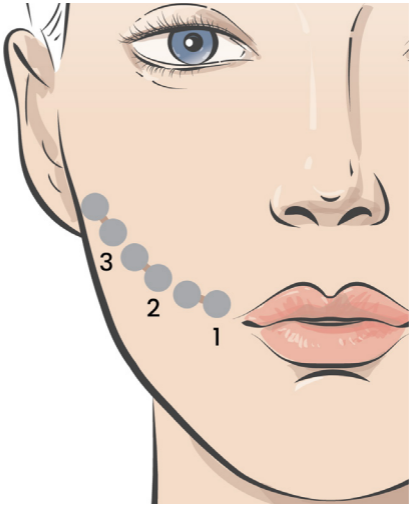- Position the Uniwand under the apple of the cheek close to the nose. Lift the cheek contour by pushing the Uniwand upwards. Hold the Uniwand in this position for 5 seconds.
- Repeat at position 2 (under the apple of the cheek, mid-way between nose and ear) and 3 (under the apple of the cheek, near the hairline).
- Repeat for the opposite side of the face.

08. Using Uniwand

1. Clean the area of the face you wish to treat.
2. Connect the Uniwand lead into either lead socket on the Uniglo device.
3. Apply the Unigel Hyaluronic Conductive Gel to the facial area. Use a similar amount as if you are applying moisturiser to your face.
4. Turn on the Uniglo device.
5. Press the programme button until FACE 2 programme has been selected. Once selected, the treatment session will automatically start. 
6. Perform the facial strokes for the selected area of the face, as indicated in the facial strokes diagrams. Each facial stroke should take 5 seconds to perform.

For each phase of the Uniglo programme, the user should complete a set of facial stroke cycles.

Phase 1 – 2 cycles

Phase 2 – 3 cycles

Phase 3 – 2 cycles

- 1 stroke = a facial stroke from one starting position
- 1 cycle = facial strokes from each starting position (see facial stroke diagrams), for both sides of the face
- 1 phase = one part of the facial treatment with a set number of cycles to complete.

For example, for the forehead, user would complete strokes on the left side of the face at position 1, 2 and 3. They would then perform the strokes on the right side of the face, again at position 1, 2 and 3. This would be one complete cycle. For phase 1, the user should complete this process twice (2 cycles), 3 times for phase 2 (3 cycles) and twice for phase 3 (2 cycles).

Note: Each phase has a preset time (see Uniglo manual) to perform the necessary cycles. However, it is expected you will complete the facial cycles for each phase within the allotted time.

7. If you have completed the number of cycles for each phase before the end of the phase time, unlock the device by pressing either ▼ button and hold the phase skip button for 3 seconds to move onto the next phase. 

Note: Expected treatment time per facial area is 5 minutes.

8. Once all three phases are complete, move onto the next area of the face you wish to treat and complete phase 1, 2 and 3 as instructed above.

Caution: Microcurrent stimulation can be uncomfortable or painful if insufficient conductive gel is used. If you start to feel pain midway through a session, apply extra conductive gel to the area of the face you are treating.

Frequency of Use

We recommend one treatment session daily, 5 days/week, for 12 weeks. You may notice some positive effects after just one session, but normally it will take

about 6 weeks before you notice a difference.

After your session

- Turn off your Uniglo device pressing and holding the power button on the side and disconnect the Uniwand.
- Wash your face with warm soapy water to remove any excess Unigel.
- It is important that the Uniwand is cleaned after each use. Clean with either an alcohol-free antibacterial wipe such as TensCare Wipes (see X-WIPES) or by wiping with a damp cloth and a solution of mild soap and water. Any excess Unigel should be removed. Wipe dry.
- If using alongside the Uniglo device, after use store the Uniwand in the Uniglo storage pouch.

09. Disposal

Clean thoroughly and dispose of according to local regulations.

10. Troubleshooting

If your Uniwand is not working properly, please check the following:

Problem	Possible Causes	Solution
Pain during treatment	Wrong programme is being used.	Uniwand should only be used with Microcurrent programmes. If using with Uniglo, this should be FACE 2.
	Not enough Unigel applied to the face.	Evenly apply a generous amount of Unigel to the areas of the face you are stimulating.
	Irritation or allergic reaction to the Unigel Hyaluronic Conductive Gel.	Try alternative water-based conductive gel. Discontinue treatment and seek advice from your healthcare professional.
No sensation	Microcurrent programme is being used.	It is normal to feel no sensation during use.

If the above review has failed to resolve your problem, or to report unexpected operation or events, or to provide feedback call TensCare or your local supplier or distributor (address on back cover) for advice.

Contact TensCare customer service on +44 (0) 1372 723 434. Our staff are trained to assist you with most issues you may have experienced, without the need to send your product in for service.

11. Warranty

Uniwand has a one-year warranty from the date of purchase. The warranty does not cover the Unigel Hyaluronic Conductive Gel.

This warranty is void if the product is modified or altered, is subject to misuse or abuse; damaged in transit; lack of responsible care; is dropped; if the unit has been immersed in water; if damage occurs by reason of failure to follow the written instruction booklet enclosed; or if product repairs are carried out without authority from TensCare Ltd.

We will repair, or at our option replace free of charge, any parts necessary to correct material or workmanship, or replace the entire unit and return to you during the period of the warranty. Otherwise, we will quote for any repair which will be carried out on acceptance of our quotation. The benefits conferred by this warranty are in addition to all other rights and remedies in respect of the product, which the consumer has under the Consumer Protection Act 1987.

Our goods come with guarantees that cannot be excluded under the UK consumer Law. You are entitled to have the goods repaired or replaced if the goods fail to be of acceptable quality.

Should repair be needed within the warranty period, please return the unit with the following information:

- Name
- Address
- Phone number
- Email address
- Proof of purchase

- Retailer's name and address
- Brief description of the problem you are experiencing

Please send the unit to:

TensCare Ltd

9 Blenheim Road

Epsom

Surrey

KT19 9BE, U.K.

Should you require any further information do not hesitate to contact us by calling +44 (0) 1372 723 434.

12. Technical Specification

Device Specifications

Dimensions 134 x 25 x 25mm

Weight 53g

Microcurrent
Parameter Range 100 - 750 μ A

Environmental Specifications

Operating Temperature range: 0-40°C
Humidity: 10-93% RH non-condensing
Atmospheric Pressure: 700-1060hPa

Storage Temperature range: -25-70°C
Humidity: 10-93% RH non-condensing
Atmospheric Pressure: 700-1060hPa

13. Symbols Table

	This symbol on the unit means "Refer to instructions for use".
	CE Mark
	Manufacturer Symbol
	Lot Number: indicates the manufacturer's batch code so that the batch or lot can be identified.
	Catalogue Number: indicates the manufacturer's catalogue number so that the device can be identified.
	Do not dispose in household waste.
	Temperature Limitation: indicates the temperature limits to which the device can be safely exposed.
	Humidity Limitation: indicates the humidity limits to which the device can be safely exposed.
	Atmospheric Pressure: indicates the atmospheric limits to which the device can be safely exposed.
	A Contraindication is used when a device should not be used because the risk of use clearly outweighs any foreseeable benefits and may result in serious injury or death.
	A Warning is used when failure to follow the instructions may result in serious injury or death.
	A Caution is used when failure to follow the instructions may result in a minor or moderate injury, or damage to the device or other property.

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TensCare Limited