Vaginal probe for electrotherapy unit

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INSTRUCTIONS FOR USE
READ CAREFULLY BEFORE USE

INDICATIONS FOR USE

The **Liberty** probe is specifically designed to suit your body. The small dimensions, smooth outer casing and curved edges mean it is easy to insert and comfortable to use.

The probe has a curved tip for easy insertion, and a flange at the base for easy removal.

For use only with Muscle Stimulation and EMG Biofeedback equipment approved to EN 60601-1 in the treatment of Stress, Urge and Mixed urinary incontinence, and pelvic pain.

For Women: The **Liberty** Probe is intended for the treatment of stress, urge and mixed urinary incontinence, and pelvic pain.

PROBE CAUTIONS

Caution: This probe is intended for single patient use only. Do not share your probe with anyone else. *Improper treatment or cross-infection may occur.*

Caution: It is important that the probe is cleaned before and after each use. *Ineffective cleaning may lead to irritation or infection.*

Caution: Never insert or remove the probe unless the control unit is powered OFF as insertion or removal when stimulation is active may cause discomfort or tissue irritation.

Caution: If tissue irritation occurs, discontinue treatment immediately. *Ask your healthcare professional for advice before continuing further treatment to prevent injury.*


Caution: Do not use a silicone-based lubricant on the metal plates of the probe as it may decrease the effectiveness of the electrical stimulation.

Caution: Electrode surface area is 8.58cm². At high output settings on some devices, this could result in a current density greater than 2mA RMS/cm². High current density may be associated with localised tissue damage or minor non-thermal burns. Increase intensity carefully and reduce if you experience sharp pain during use or discharge or bleeding after use.

Warning: Do not use if the pack is opened or damaged.

INSTRUCTIONS FOR USE

1. Clean the probe with an alcohol-free anti-bacterial wipe or with warm soapy water.
2. Connect the lead from the base of the unit to the lead in the probe. Push the pin ends firmly into the pigtail ends of the probe lead.
3. Before using the probe, you will need to visit the toilet.
4. Lubricate the metal electrode surfaces and probe tip with a water-based lubricant, such as TensCare Go Gel or water.

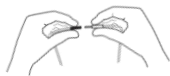-  **Caution:** Do not use a silicone based lubricant on the stimulation contacts as it may decrease the effectiveness of the muscle stimulation.
5. Choose a comfortable position, such as lying down on your bed on your side with your knees raised.

-  **Warning:** Ensure the unit is switched OFF before inserting the probe.

6. **Vaginal insertion:** After the wire is securely connected, insert the probe into your vagina, in the same way as you would a tampon, until only the plastic ring at the base of the probe is visible. The metal parts conduct the electrical pulse and should be in contact with the main part of the muscle at all times. The tissues close to the entrance are more sensitive, so you should avoid stimulating them. Probes with long electrodes (the metal part) that run up and down the length of the attachment should always be inserted with the metal parts facing hip-to-hip.
7. After your training session, switch the control unit off, then remove the probe by holding the positioning end rim and gently pulling outwards.

CLEANING

It is important that the probe is cleaned before and after each use. Clean with either an alcohol-free antibacterial wipe such as TensCare Wipes (see X-WIPES) or by wiping with warm soapy water. Rinse under running water and dry thoroughly and return the unit to the storage pouch.

Warnings:

Do not immerse the probe in a liquid.

Keep away from children under 5 years of age. Long cable - risk of strangulation in infants.

Keep the device away from lint and dust, as long-term exposure to lint or dust may affect the socket.

Keep the device away from sunlight, as long-term exposure to sunlight may affect the rubber causing it to become less elastic and crack.

DISPOSAL

Clean thoroughly and dispose of according to local regulations.


TECHNICAL SPECIFICATIONS

Material: Electrode body: Plastic
ABS PA-757

Electrodes: Stainless steel #304

Cables: PVC over copper tinsel
conductor

“In Use” product lifetime: 6
months

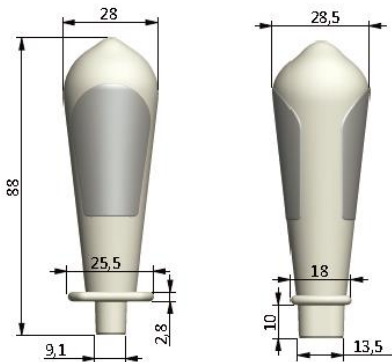ENVIRONMENTAL CONDITIONS:

Storage: +5 to +40°C, RH up to 93%, 700 to 1060 hPa

Operation: +5 to +40°C, RH 15 to 93%, 700 to 1060 hPa

Water ingress protection: IP22

The first number 2: Protected against access to hazardous parts with a finger, and the jointed test finger of 12 mm \varnothing , 80 mm length

The second number 2: Protected against vertically falling water drops when enclosure is tilted up to 15°.

This symbol on the device means “Refer to instructions for use”

To report unexpected operation or events contact Tenscare Ltd.

Please find our details on the final page of this insert.

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TensCare Ltd,
Advena Ltd.